

Hella honeydew

Sample ID: BIA260210S0178
Strain: MANU0041-0359-1-HH-26-1
Harvest Lot:
Matrix: Concentrates & Extracts
Type: Distillate
Sample Size: 1 units
Lot#:

Produced:
Collected:
Received: 02/10/2026
Completed: 02/19/2026
Batch#:

Client:
Blondie's Bud company
Lic. # CLTV0359
 232 Neshobe Cir
 Brandon, VT 05733



Summary

| Test | Date Tested | Result |
|--------------|-------------|----------|
| Sample | | Complete |
| Cannabinoids | 02/12/2026 | Complete |
| Terpenes | 02/10/2026 | Complete |

Cannabinoids

Completed

| 84.97% Total THC | | | | | | 0.35% Total CBD | | 91.71% Total Cannabinoids | | | | | |
|---------------------|--------|---------|---------|-------|--------------|--------------------|--------|------------------------------|---------------|-------------|--------------|--|--|
| Analyte | LOQ | Results | Results | Mass | Mass | Analyte | LOQ | Results | Results | Mass | Mass | | |
| | mg/g | % | mg/g | mg/mL | mg/container | | mg/g | % | mg/g | mg/mL | mg/container | | |
| CBDVa | 0.0003 | <LOQ | <LOQ | | | CBCVa | 0.0003 | <LOQ | <LOQ | | | | |
| CBDV | 0.0003 | <LOQ | <LOQ | | | CBNa | 0.0003 | <LOQ | <LOQ | | | | |
| CBDa | 0.0005 | <LOQ | <LOQ | | | Δ9-THC | 0.0005 | 84.97 | 849.7 | | | | |
| CBGa | 0.0005 | <LOQ | <LOQ | | | Δ8-THC | 0.0003 | <LOQ | <LOQ | | | | |
| CBG | 0.0005 | 3.48 | 34.8 | | | Δ10-THC* | 0.0002 | <LOQ | <LOQ | | | | |
| CBD | 0.0005 | 0.35 | 3.5 | | | CBL | 0.0005 | <LOQ | <LOQ | | | | |
| THCV | 0.0003 | 0.72 | 7.2 | | | CBC | 0.0003 | 1.07 | 10.7 | | | | |
| CBLV | 0.0003 | <LOQ | <LOQ | | | THCa | 0.0005 | <LOQ | <LOQ | | | | |
| CBCV | 0.0003 | <LOQ | <LOQ | | | CBCa | 0.0006 | <LOQ | <LOQ | | | | |
| THCVa | 0.0003 | <LOQ | <LOQ | | | CBLa | 0.0005 | <LOQ | <LOQ | | | | |
| CBN | 0.0005 | 1.12 | 11.2 | | | Total THC | | 84.97 | 849.75 | | | | |
| | | | | | | Total CBD | | 0.35 | 3.46 | | | | |
| | | | | | | Total | | 91.71 | 917.09 | 0.00 | 0.00 | | |

Analyst: 048

Cannabinoids Methodology: High Performance Liquid Chromatography (HPLC) using PerkinElmer FLEXAR™ with Photo Diode Array Detector (PDA)

Total CBD and total THC are calculated values, to account for assumed decarboxylation from the acid form (THCA or CBDA) to the neutral form, causing weight loss of the acid group. These values are calculated as follows:

$$\text{Total THC} = (\text{THCA} \times 0.877) + \Delta 9\text{-THC}$$

$$\text{Total CBD} = (\text{CBDA} \times 0.877) + \text{CBD Reagent}$$

Blanks: < LOQs for all analytes

LOQ = The lowest quantity that this method can reliably detect. Any cannabinoid that was not detected is assumed to be less than the stated LOQ (<LOQ).

All results reflect dry weight of material, based on % moisture of the sample.

Measurement of Uncertainty (MU): the parameter, associated with the result of a measurement, that characterizes the dispersion of the values that could reasonably be attributed to the

particular quantity subject to measurement. Δ9-THC MU = ±0.005% Total THC MU = ±0.007%

All other cannabinoid MU values are available upon request.

All moisture and water activity analysis is determined by dewpoint measurement using an AQUALAB water activity meter.

*The result is the sum of delta-10 isomers.




Luke Emerson-Mason
 Laboratory Director
 02/19/2026

Confident LIMS
 All Rights Reserved
coa.support@confidentlims.com
 (866) 506-5866
www.confidentlims.com



Hella honeydew

Sample ID: BIA260210S0178
Strain: MANU0041-0359-1-HH-26-1
Harvest Lot:
Matrix: Concentrates & Extracts
Type: Distillate
Sample Size: 1 units
Lot#:

Produced:
Collected:
Received: 02/10/2026
Completed: 02/19/2026
Batch#:


Client:
Blondie's Bud company
Lic. # CLTV0359
 232 Neshobe Cir
 Brandon, VT 05733

Terpenes

Completed

| Analyte | LOQ | Results | Results |
|---------------------|-------|---------------|--------------|
| | mg/g | mg/g | % |
| Ocimene | 0.010 | 4.899 | 0.490 |
| Limonene | 0.010 | 4.152 | 0.415 |
| β-Myrcene | 0.010 | 4.039 | 0.404 |
| α-Pinene | 0.010 | 4.012 | 0.401 |
| β-Pinene | 0.010 | 3.420 | 0.342 |
| β-Caryophyllene | 0.010 | 3.222 | 0.322 |
| Linalool | 0.010 | 2.140 | 0.214 |
| Terpinolene | 0.010 | 1.614 | 0.161 |
| α-Humulene | 0.010 | 1.033 | 0.103 |
| Camphene | 0.010 | 0.371 | 0.037 |
| Eucalyptol | 0.010 | 0.174 | 0.017 |
| 3-Carene | 0.010 | 0.168 | 0.017 |
| γ-Terpinene | 0.010 | 0.127 | 0.013 |
| α-Terpinene | 0.010 | 0.078 | 0.008 |
| α-Bisabolol | 0.010 | 0.028 | 0.003 |
| Caryophyllene Oxide | 0.010 | 0.027 | 0.003 |
| cis-Nerolidol | 0.010 | <LOQ | <LOQ |
| Geraniol | 0.010 | <LOQ | <LOQ |
| Guaiol | 0.010 | <LOQ | <LOQ |
| Isopulegol | 0.010 | <LOQ | <LOQ |
| p-Cymene | 0.010 | <LOQ | <LOQ |
| trans-Nerolidol | 0.010 | <LOQ | <LOQ |
| Total | | 29.504 | 2.950 |

Primary Aromas

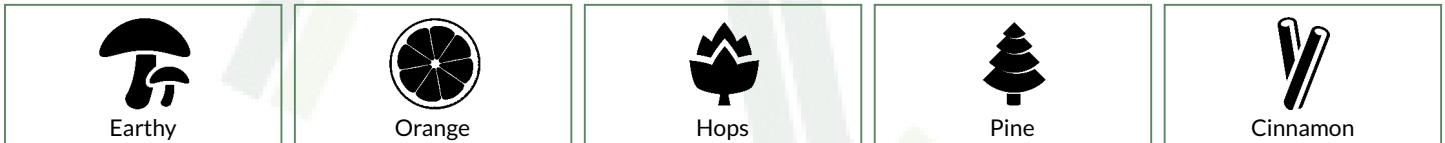

Analyst: 063

LOQ = The lowest quantity this method can reliably detect. Any terpene that was not detected is assumed to be less than the stated LOQ (<LOQ).

Terpene Methodology: Headspace Sampler, Gas Chromatography-Mass Spectrometry (GC-MS), using Perkin Elmer Clarus® SQ8 GC MS

Reagent Blanks: < LOQs for all analytes

All results reflect dry weight of material, based on % moisture of the sample.

All moisture and water activity analysis is determined by dewpoint measurement using an AQUALAB water activity meter.




Luke Emerson-Mason
 Laboratory Director
 02/19/2026

Confident LIMS
 All Rights Reserved
coa.support@confidentlims.com
 (866) 506-5866
www.confidentlims.com

